



YOUR VISION. OUR FOCUS.

# Enhanced Auto Kerato-Refractometer with Corneal Mapping


## KR-800PA

Auto Kerato-Refractometer



### Corneal Mapping

The Corneal mapping function provides a patient's corneal curvature radius and refraction power map. By performing corneal mapping in conjunction with the refraction measurement, it is possible to observe abnormal corneal shape with total refraction.

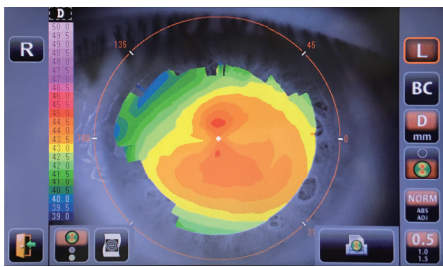


### Unparalleled Accuracy

The unique\* Rotary Prism Technology, exclusive to TOPCON, allows for unparalleled precision and reliability. This system measures a wider area of the retina while permitting measurements through pupils as small as 2 mm in diameter.

\*The off-centered ring target can cover areas that may be occluded by a small pupil. The angle of the prism enables a wider area of the retina to be measured.

### Easy Operated 3D Auto Alignment Function

The 3D auto alignment feature supports ease of use by maintaining XYZ alignment even during small eye movements, providing fast and precise alignment. An operator follows the alignment guidance on the monitor for initial alignment, then 3D auto alignment and measurement start automatically.

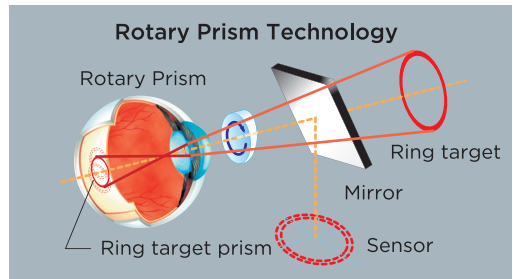

※image

\* Subject to change in design and/or specifications without advanced notice.

**IMPORTANT**

In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.



**TOPCON CORPORATION**

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan. Phone:+81-(0)3-3558-2522/2502 Fax:+81-(0)3-3965-6898 www.topcon.co.jp

#### TOPCON MEDICAL SYSTEMS, INC.

111 Bauer Drive, Oakland, NJ 07436, U.S.A.  
Phone:+1-201-599-9100 Fax:+1-201-599-5250 www.topconmedical.com

#### TOPCON HEALTHCARE SOLUTIONS, INC.

111 Bauer Drive, Oakland, NJ 07436, U.S.A.  
Phone:+1-201-599-5100 Fax:+1-201-599-5250 www.topconhealth.com

#### TOPCON CANADA INC.

110 Frovicher Avenue, Boisbriand, QC J7G 1N1 CANADA  
Phone:+1-450-430-7771 Fax:+1-450-430-6457 www.topcon.ca

#### TOPCON EUROPE MEDICAL B.V.

Essebaan 11, 2908 LJ Capelle a/d IJssel, P.O.Box145,  
2900 AC Capelle a/d IJssel, THE NETHERLANDS  
Phone:+31-(0)10-4585077 Fax:+31-(0)10-4585045  
E-mail: medical@topcon.eu; www.topcon-medical.eu

#### ITALY OFFICE

Viale dell'Industria 50, 20037 Paderno Dugnano, (Milano), ITALY  
Phone:+39-02-91866071 Fax:+39-02-91081091 E-mail: info@topcon.it; www.topcon.it

#### DANMARK OFFICE

Praestemarksvej 25, 4000 Roskilde, DANMARK  
Phone:+45-46-327500 Fax:+45-46-327555  
E-mail: topcon@topcondanmark.dk www.topcondanmark.dk

#### IRELAND OFFICE

Unit 276, Blanchardstown, Corporate Park 2 Ballycoolin Dublin 15, IRELAND  
Phone:+353-18975900 Fax:+353-18233915 E-mail: medical@topcon.ie; www.topcon.ie

#### TOPCON FRANCE MEDICAL S.A.S.

BAT AI, 3 route de la revolte, 93206 Saint Denis Cedex, France  
Phone:+33-1-49212323 Fax:+33-1-49210494 E-mail:topcon@topcon.fr; www.topcon.fr

#### TOPCON DEUTSCHLAND MEDICAL G.m.b.H.

Hanns-Martin-Schleyer-Strasse 41, D-47877 Willich, GERMANY  
Phone:+49-(0)2154-8850 Fax:+49-(0)2154-88517 E-mail:med@topcon.de; www.topcon.de

#### TOPCON SCANDINAVIA A.B.

Neongatan 2, P.O.Box 25, 43151 Malmdal, SWEDEN  
Phone:+46-(0)31-7109200 Fax:+46-(0)31-7109249 E-mail:medical@topcon.se; www.topcon.se

#### TOPCON ESPAÑA S.A.

HEAD OFFICE  
Frederic Mompou 4 Esc. A Bajos 3, 08960 Sant Just Desvern Barcelona, SPAIN  
Phone:+34-93-4734057 Fax:+34-93-4733932 E-mail: medica@topcon.es; www.topcon.es

#### TOPCON (GREAT BRITAIN) MEDICAL LTD.

Topcon House, Kennet Side, Bone Lane, Newbury, RG14 5PX U.K.  
Phone:+44-1635-551120 Fax:+44-1635-551170 E-mail:medical@topcon.co.uk; www.topcon.co.uk

#### TOPCON POLSKA Sp. z o.o.

ul. Warszawska 23, 42-470 Siewierz, POLAND  
Phone:+48-(0)32-6705045 Fax:+48-(0)32-6713405 www.topcon-polska.pl

#### Mehra Eyetech Private Limited

801 B Wing, Lotus Corporate Park, Graham Firth Steel Compound Goregaon (East)  
Mumbai 400063 Maharashtra, INDIA  
Phone:+91-22-61285455 www.mehraeyetech.in

#### TOPCON SINGAPORE MEDICAL PTE. LTD.

1 JALAN KILANG TIMOR #09-01 PACIFIC TECH CENTRE SINGAPORE 159303  
Phone:+65-68720806 Fax:+65-67735550  
E-mail:medical\_sales@topcon.com.sg www.topcon.com.sg

#### REPRESENTATIVE OFFICE IN INDONESIA

Level 38, Tower A, Kota Kasablanka unit GH-04, Jl.  
Casablanca, Kav 88, Jakarta, Indonesia 12870  
Phone:+62-21-2963-8004

#### TOPCON INSTRUMENTS (MALAYSIA) SDN.BHD.

No. D1, (Ground Floor), Jalan Excella 2, Off Jalan Ampang Putra,  
Taman Ampang Hilir, 55100 Kuala Lumpur, MALAYSIA  
Phone:+60-(0)3-42709866 Fax:+60-(0)3-42709766

#### TOPCON INSTRUMENTS (THAILAND) CO.,LTD.

77/162 Simsathorn Tower, 37th Floor, Krungthoburi Rd., Klongtonnai,  
Klongsarn, Bangkok 10600, THAILAND  
Phone:+66-(0)2-440-1152-7 Fax:+66-(0)2-440-1158

#### TOPCON (BEIJING) MEDICAL TECHNOLOGY CO., LTD.

Room 2808, Tower C, JinChangAn Building No.82 Middle Section of East 4th Ring Road,  
Chaoyang District, Beijing 100124, P.R.China  
Phone:+86-10-87945176



TOPCON CORPORATION



EC9711081  
TOPCON CORPORATION



YOUR VISION. OUR FOCUS.



new

# KR-800PA

Auto Kerato-Refractometer

*Enhanced Auto Kerato-Refractometer  
with Corneal Mapping*

- » Corneal Mapping Function
- » Topcon's Rotary Prism Technology for Unparalleled Accuracy
- » Easy Operation with 3D Auto Alignment
- » Compact Design
- » Wide Touchscreen Control Monitor
- » Simple and Intuitive Data Display

